

Invacare® | **Matrix® Libra®**

NEW



The next generation cushion - Superior comfort with the highest level of protection and positioning

- Improved design combines foam and fluid to provide the highest level of skin protection, positioning and adjustability for unique user needs
- New dual layer fluid sac using Invacare's **Flo-tech™** technology for enhanced shear management and pressure redistribution
- Specially contoured to provide balance and stability to the pelvis, optimizing body function
- Lightweight and maintenance free, the **Libra** balances both functional needs with an active lifestyle
- Recommended for all users including those with specific individual seating needs

matrix®
SEATING SERIES



Yes, you can.®



Unique design

- Foam base – anatomically contoured to optimize pelvic support and positioning
- Durable high resilient foam made with Ultra-Fresh™ anti-microbial and odor protection
- Fluid sac – Dual layer with sealed compartments for maximum protection around bony prominences

Functional benefits

- Pre-ischial support prevents sliding to reduce shear
- Trochanteric shelf promotes level pelvis for lateral stability
- Extended rear cushion radius helps prevent posterior pelvic tilt
- Waterfall front edge accommodates tight hamstrings and foot propulsion
- Ab/adductor contour for comfortable midline leg position

Reliable and lightweight

- Self-equalizing fluid is non-temperature sensitive. No kneading or maintenance is required
- With a weight reduction of 40% compared to other high-end **Flo-tech™** cushions, the **Matrx Libra** is lightweight and suitable for an active lifestyle
- For complex seating needs, a range of easy to use positioning accessories are available.

Features and Options



Envelopment

Increased surface area for superior pressure redistribution.



Pressure Re-distribution

Load transferred from pelvis to thighs reduces pressure in at-risk areas.



Immersion

Bony prominences comfortably immerse into fluid sac for optimal protection.



Shear Reduction


Unique design allows controlled pelvis movement without migration.

Features and Options



Waterfall cut

For tight hamstrings, easier transfers and foot propulsion.

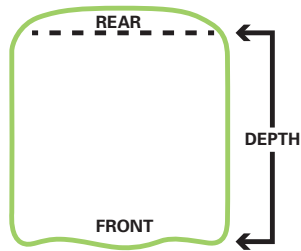

Cushion size

Seat depth is measured along the side of the cushion.



Rear cushion radius

A 3.8 cm curved extension for maximum surface area minimizes gap between seat and backrest. Stabilizes fluid in tilted positions, with coccyx cut out for added relief.



Anterior wedge

Increases slope, prevents sliding and reduces shear.



Trimmed anterior wedge

For hip flexion limitations and/or lower seat height for foot propulsion.



Tapered wedge

For enhanced adductor support.



Tapered wedge

For obliquity and pelvic rotation when used in rear.



Lateral wedge

Reinforces lateral stability for enhanced midline positioning.



Fluid supplement sac

For targeted protection/obliquity accommodation. Fits securely inside main sac and pelvic well.



Abductor

For comfortable mid-line positioning of thighs and hips.



Side pocket for convenient storage

Includes quad pull cord.

Available sizes

Standard version in cm: width x depth Maximum user weight: 136 kg

36 x 31	38 x 41	41 x 41	43 x 41	46 x 41	48 x 41	51 x 41
36 x 36	38 x 46	41 x 46	43 x 46	46 x 46	48 x 46	51 x 46
36 x 41	38 x 51	41 x 51	43 x 51	46 x 51	48 x 51	51 x 51

Heavy duty version in cm: width x depth Maximum user weight: 227 kg

53 x 46	56 x 46
53 x 51	56 x 51
53 x 56	56 x 56

Technical data

Matrx Libra



Weight: 1,6 kg



Max User Weight:
Std 136 kg
HD 227kg



Outer cover:
machine wash
80°




Inner liner:
surface clean
with damp cloth




- 1/ Compliant with the directive on the Restriction of the use of certain Hazardous Substances in electrical and electronic equipment – RoHS (2011/65/EC).
- 2/ Compliant with the regulation concerning the Registration, Evaluation, Authorization and Restriction of Chemicals.

Cover options

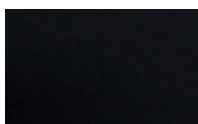

Black Smooth
Water resistant
Vapour permeable

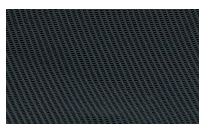

Spacer

Invacare International Sàrl
Route de Cité-Ouest 2 - 1196 Gland
Switzerland
Tel: +41 (0)22 354 60 10
Fax: +41 (0)22 354 60 11
hqeurope@invacare.com
www.invacare.eu.com

© 2015 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print.
Invacare reserves the right to alter product specifications without prior consultation.

Matrx Libra - EU - 02/2015



Yes, you can.®